

DS LABORATORY ROTATOR/ STIRRER

DSR SERIES


INSTRUCTION MANUAL:

Operation is very easy to do....

SPEED CONTROLLER

*Turn to desired setting for slow or fast shaking or stirring]

*always start at slow or lowest number setting. Gradually increasing until desired action is attained

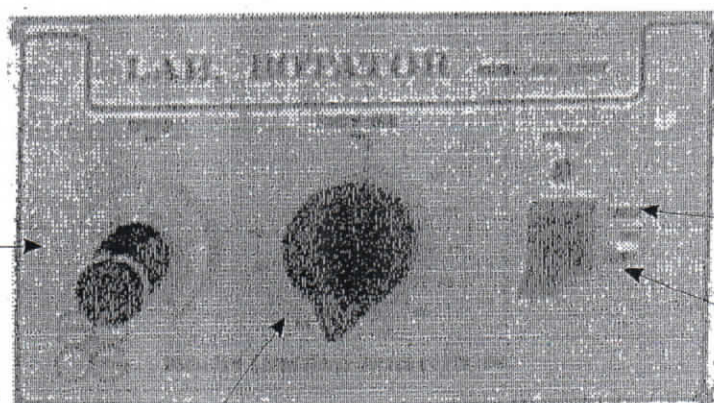

OPERATION SELECTOR


TIMED OPERATION
[FLICK UPWARD]

CONTINUOUS OPERATION
[FLICK DOWNWARD]

TIMER [PLEASE TURN FULLY THEN SET TO DESIRED TIME.]

* For fast rotating shaking or stirring,
Fasten or restrain container with
spring retainer

* For slower action you can eliminate
the use of retainer.

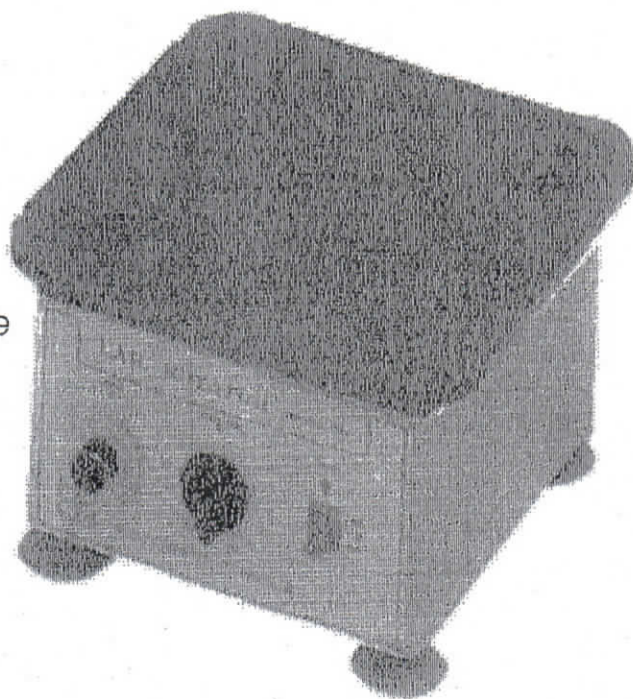

DSR-2100V

DIGISYSTEM LABORATORY INSTRUMENTS, INC.

P.O.BOX 12236 TAIPEI, TAIWAN, R.O.C.